


# DACLON

## NYLON

DACLON  
NYLON

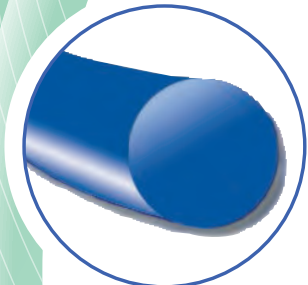

|                        |   |
|------------------------|---|
| <b>Type</b>            | Monofilament  |
| <b>Composition</b>     | Extrusion of polyamide 6.0 or 6.6   |
| <b>Coating</b>         | None  |
| <b>Colour</b>          | Blue<br>Black   |
| <b>Tissue reaction</b> | Minimal   |
| <b>Absorption</b>      | Non absorbable, gradually encapsulated by connective tissue. The thread mass diminishes, approximately 10% a year by rupture of chemical links (hydrolytic action). |
| <b>Presentation</b>    | Box with 12 sutures   |

*For quality  
&  
safety*

**smi**

[www.sutures.be](http://www.sutures.be)

# DACLON

## NYLON

### **CHARACTERISTICS**

- ☑ High tensile strength, non absorbable
- ☑ Synthetic
- ☑ Monofilament with microscopic uniform diameter
- ☑ Pliability and softness
- ☑ Hermetically sealed packing

### **BENEFITS**

- Excellent and permanent tissue support
- Minimal tissue reaction
- Soft passage through the tissues without sawing, tissue drag and trauma
- Absence of capillarity (no support of bacterial growth)
- Excellent knot tie down
- Good handling properties
- Easy and secure knotting
- Guaranteed seal and product sterility

**INDICATIONS**    General soft tissue closing and/or ligation

**Daclon Blue**  
 especially in    Skin closure  
                       General surgery  
                       Plastic surgery  
                       Gastrointestinal surgery  
                       Gynaecology  
                       Obstetrics  
                       Orthopaedics

**Daclon Black**  
 especially in    Plastic surgery  
                       Ophthalmic surgery